

Zoé Sleep App — Setup Guide

Longevity Retreat, Half Moon Bay | January 16-18, 2025

A Note Before You Begin

This is not a passive sleep tracker. This is your **personal hybrid digital-human sleep clinic experience**.

If it feels like work — good. It is. You're working on your sleep.

I aim to be scientifically accurate. The validated clinical instruments we're using require this level of engagement. There are no shortcuts. The burden is the point.

Yes, it's hard work. But it's only 10 days. After that, it gets easier. We move from assessment into the prescription phase — where you'll receive your personalized sleep protocol and the app shifts to supporting your improvement journey.

At the end, you won't get a generic report — you'll get my personal analysis and a protocol designed specifically for you.

Ten days. Ten minutes a day. Your sleep is worth it.

Quick Setup (5 Minutes)


1. Install **TestFlight** from the App Store: apps.apple.com/app/testflight
 2. Open the Zoé Sleep invite: testflight.apple.com/join/62aRAuVH
 3. Install Zoé Sleep, create account, complete onboarding
 4. Follow the 6-step guided tour
 5. Tap your **initial (top-right)** → **Health Data Settings** → **Sync Health Data**
 6. **Apple Watch users:** Open app on both devices → tap **Sync Now**
- Done.** Complete Sleep Log + Assessment daily for 10 days.

Step 1: Install TestFlight & Zoé Sleep

Download **TestFlight** from the App Store: apps.apple.com/app/testflight


Then open the Zoé Sleep invite link: testflight.apple.com/join/62aRAuVH

Tap **Accept**, then **Install**.

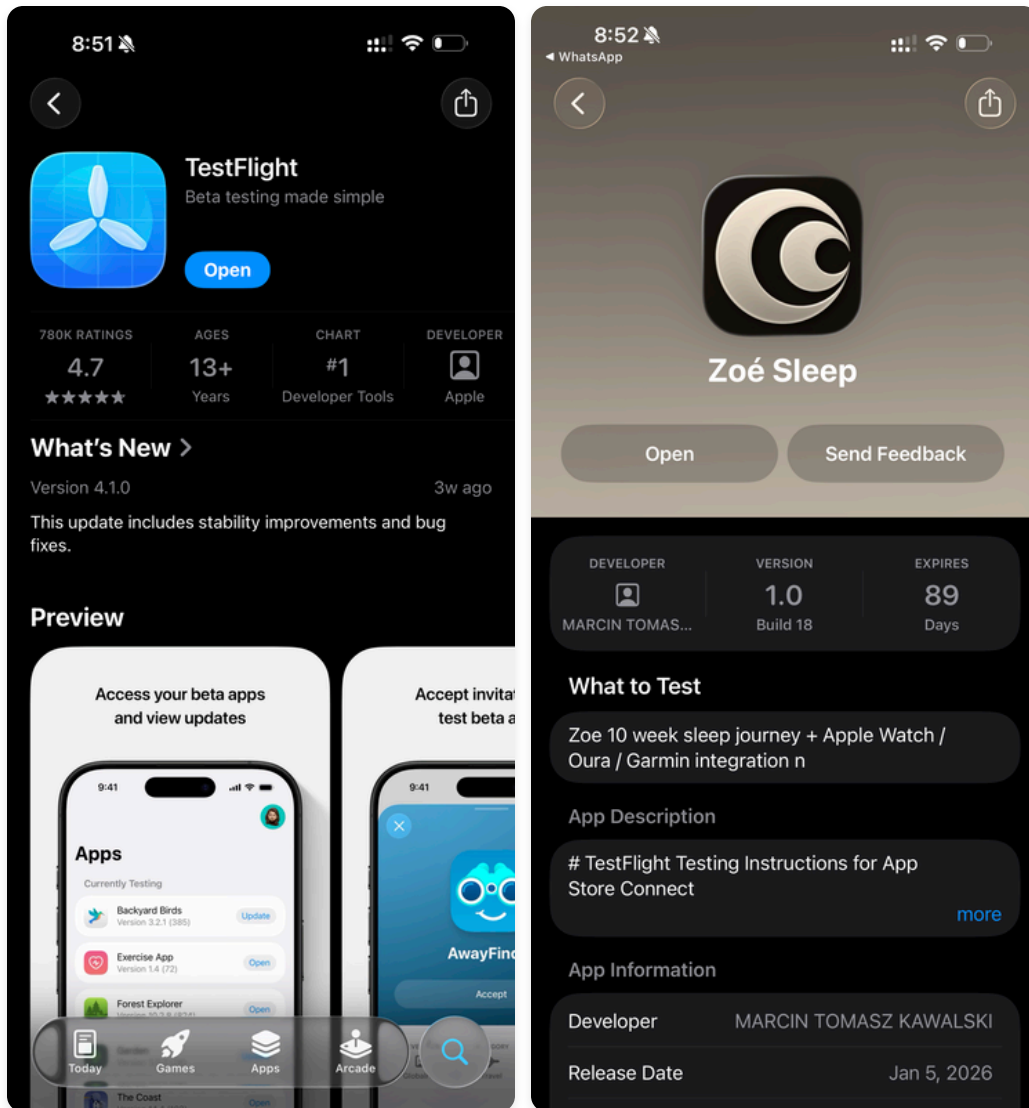


Step 2: Create Account & Onboarding

Sign up with email or Apple ID:

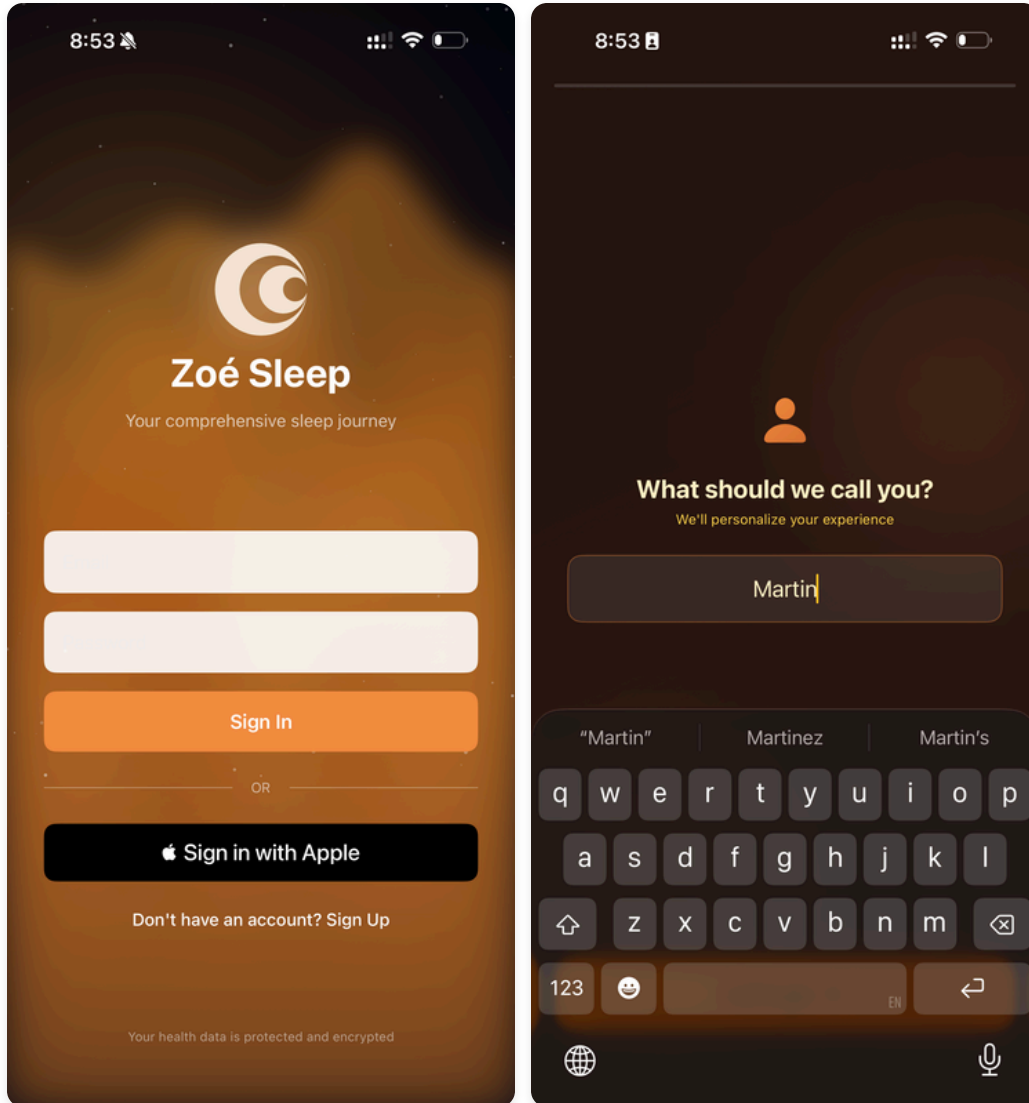


The app will guide you through onboarding — it will ask about your **height, weight, gender, and age**:

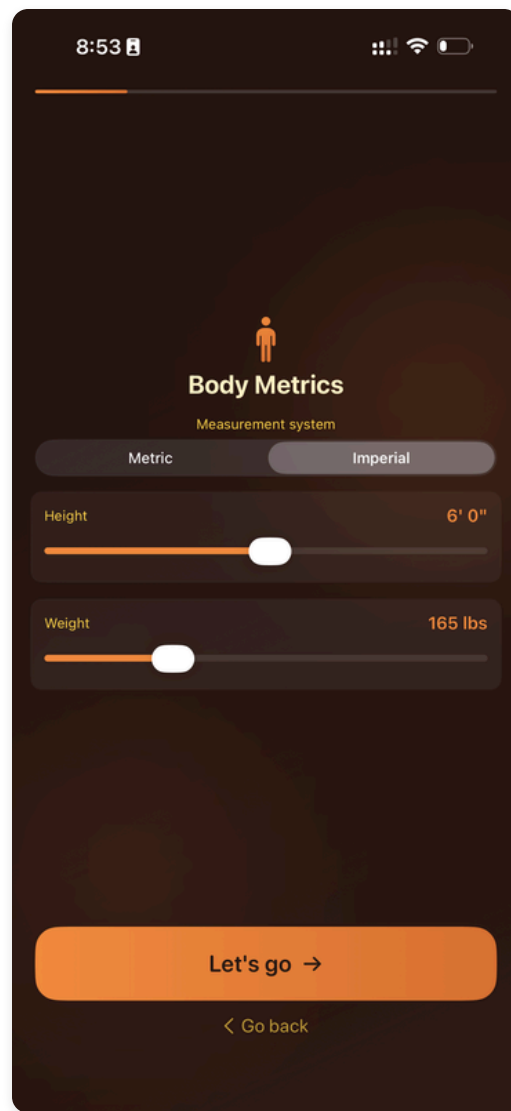

Allow HealthKit access when prompted (this grants permission — the actual data sync happens later).

Step 3: Guided Tour


After onboarding, you'll land on the main dashboard. The app automatically enters **tutorial mode** — a 6-step guided tour explaining each screen element.

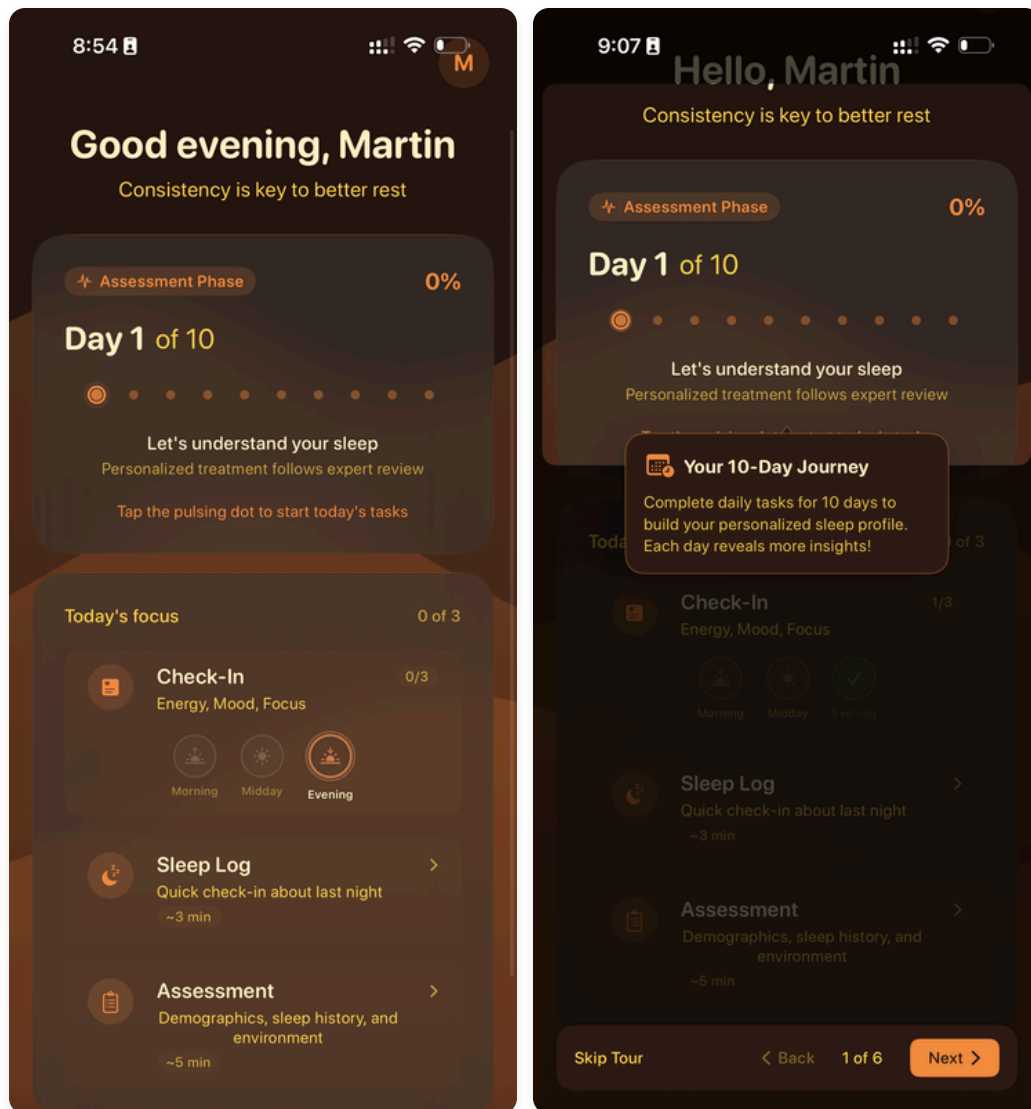

Follow along and tap **Next** to proceed through all 6 steps. The tour explains your daily tasks: Check-In, Sleep Log, and Assessment.

Step 4: Sync Your Health Data — Critical!

After the tour, tap your **initial** (top-right) → **Health Data Settings** → **Sync Health Data**

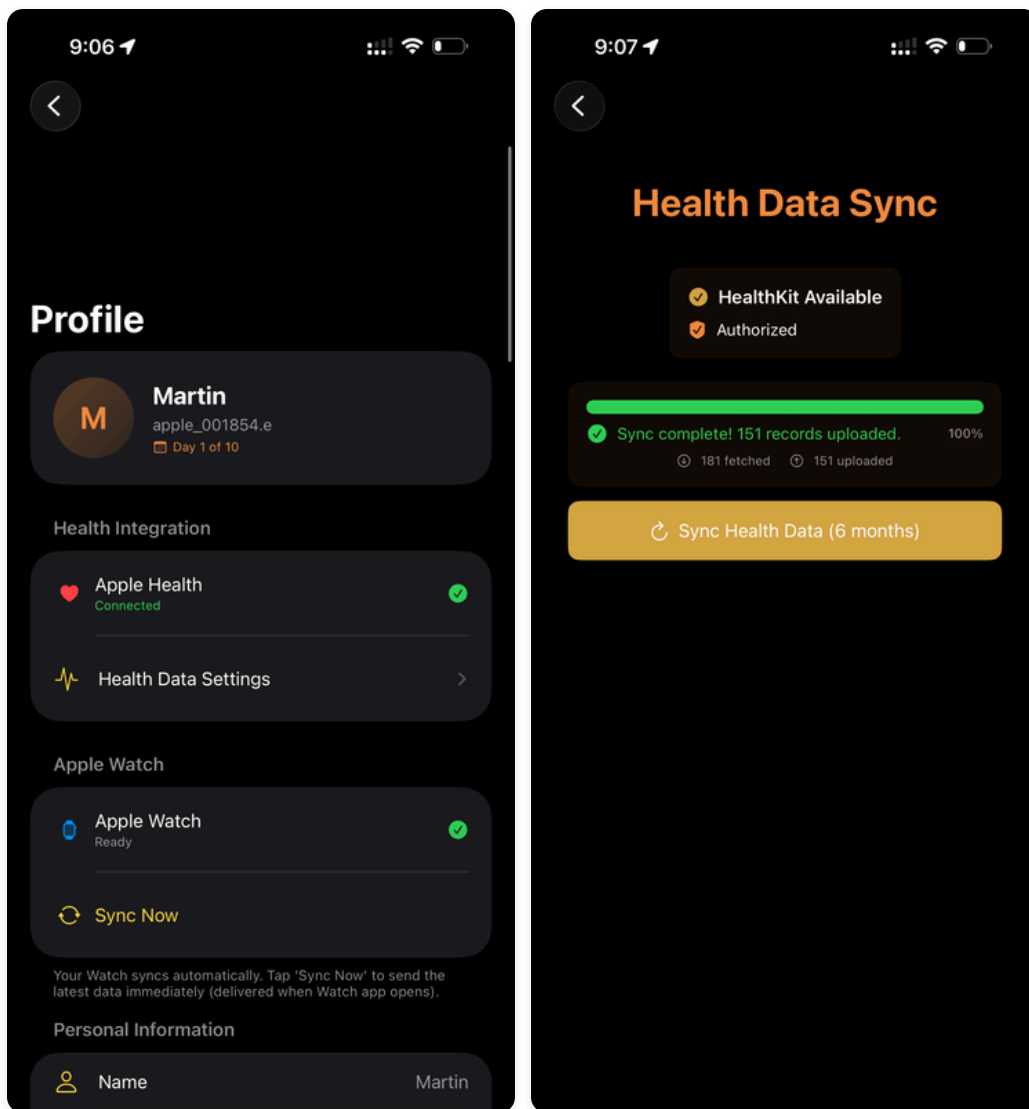


This imports data from Apple Watch, Oura, Garmin, Whoop — any device syncing to Apple Health.

Apple Watch Users Only

The Zoé Sleep Watch app installs automatically when you install the iPhone app.

To sync iPhone ↔ Watch: Both apps must be open at the same time. Then in iPhone Settings, tap **Sync Now**.

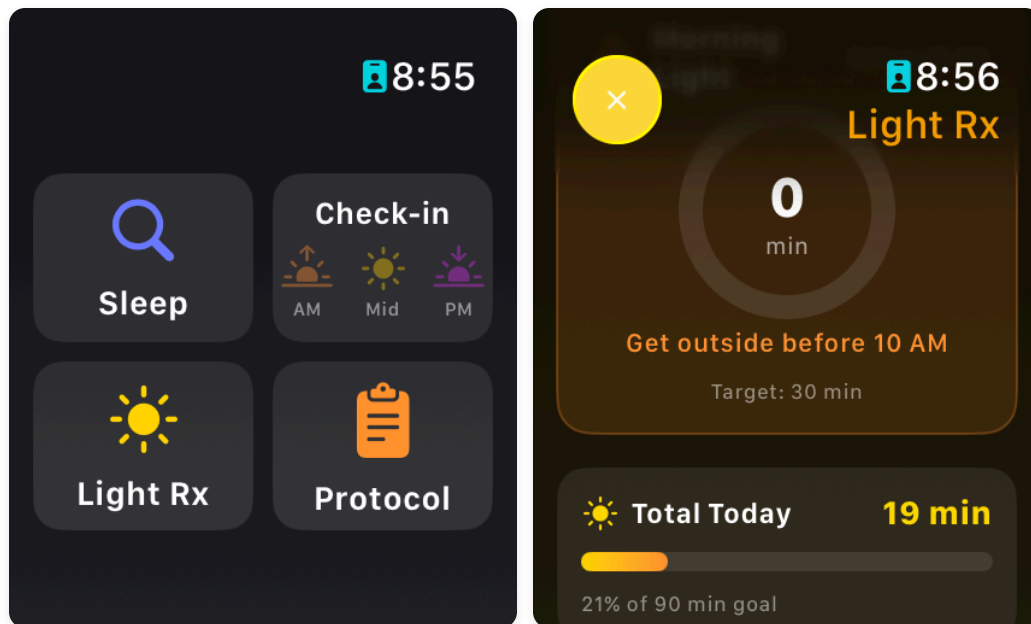

Watch features:

- **Check-in:** Use the Digital Crown to quickly rate energy/mood/focus
- **Light Rx:** Track morning light exposure (aim for 30 min before 10 AM)


Light exposure is the strongest driver of circadian rhythm. Check your progress throughout the day.

Your Daily Routine

For 10 days, complete these each morning:


Task	Time	Required
Sleep Log	~3 min	Yes
Assessment	~5 min	Yes
Check-ins (3x)	30 sec	Optional

Sleep Log

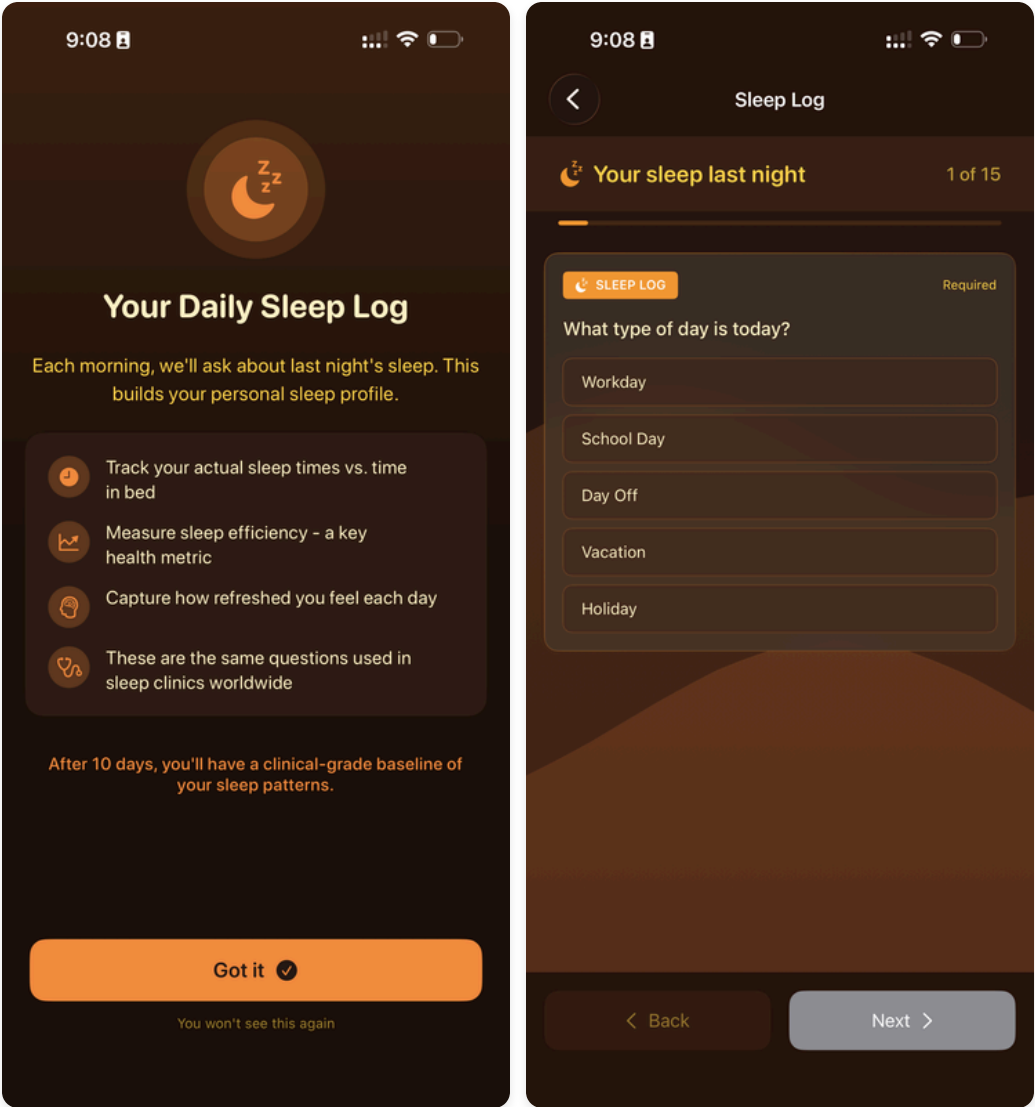

Based on the Stanford Sleep Log. Answer from **subjective memory** — don't check your wearable first! The gap between perception and reality is clinically valuable.

Assessment

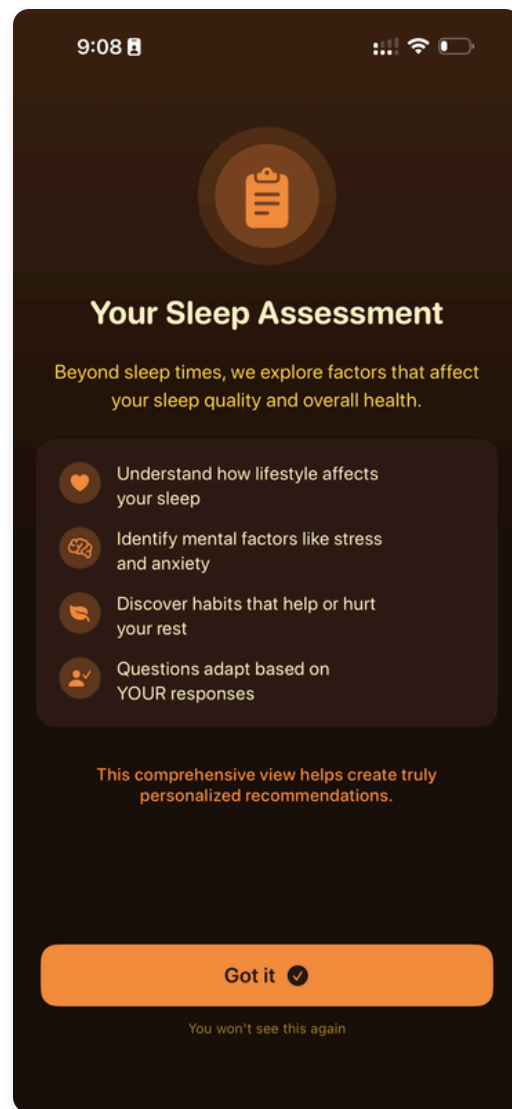


Days 1-5: Core questions (similar for everyone) **Days 6-10:** Personalized deep dives based on your answers (insomnia, sleep apnea, shift work, etc.)

Daily Lockout

App locks after daily tasks complete, unlocks at **4:00 AM** next day. This ensures data quality.

Settings Menu

Tap your initial (top-right) to access settings:

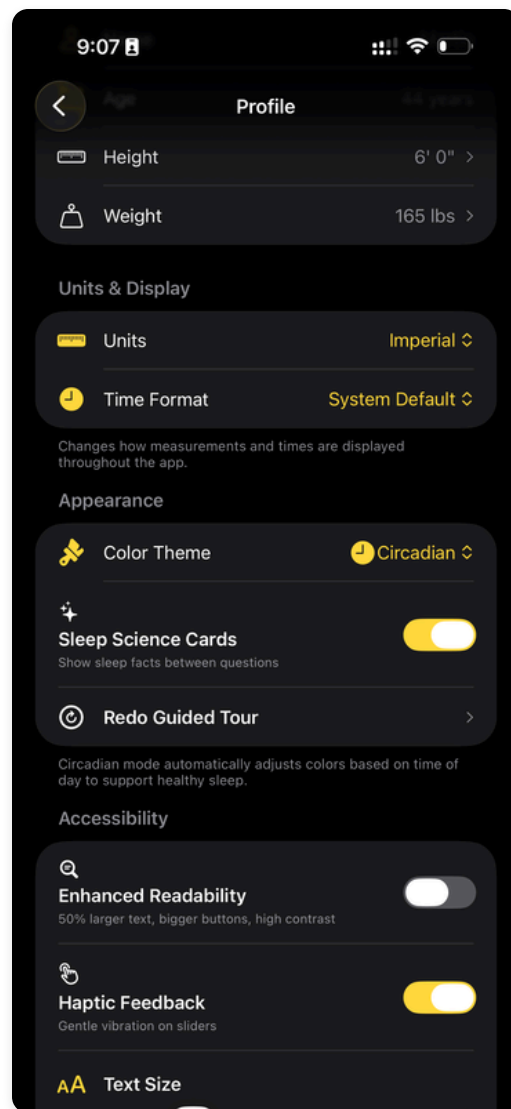

Setting

Description

Health Data Settings

Sync wearable data

Personal Information

Review/edit profile

Color Theme

Circadian mode adjusts by time of day

Sleep Science Cards

Toggle educational cards during assessments

Debug Mode

Experimental features (not fully tested)

After 10 Days

I'll personally analyze your sleep profile and questionnaire responses. At the retreat, we'll review your results together and I'll assign a personalized sleep improvement protocol.

Tips: Be consistent (same time each morning) · Don't peek at wearable data before Sleep Log · Don't skip days — bad nights are valuable · Get morning light

Questions? Contact Me

- **Email:** kawalski@stanford.edu
- **iMessage:** +1 (650) 705-9311
- **WhatsApp:** +41 78 753 3505

See you in Half Moon Bay.